

Operating datasheet

PHOTOPOD:

numerical photometer



I-010-CYANURIC ACID: 10 - 200 mg/l

Photopod

LS

Reagent kit: 1MT130
Preparation time: ~ 5min

REAGENTS

Cyanuric Acid Reagent 1RA020

EQUIPMENT

Graduated Plastic Tube 14TP00
Glass Tube 1CR099
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 30 drops of cyanuric acid reagent
Close and shake vigorously.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **010 Cya.Ac. : 10 - 200mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer.
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT301.

I-011- CYANURIC ACID: 10 - 200 mg/l

Reagent kit: 1MT048
Preparation time: ~ 5min

Photopod

SP*REAGENTS*

Cyanuric Acid Tablets 1CA008

EQUIPMENT

Graduated Plastic Tube 14TP00
Crushing Rod 1AP018
Glass Tube 1CR099
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Cyanuric Acid Tablet then wait 2 minutes dissolution
Close the tube and shake vigorously for 2 min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 1 minute.
Proceed to the measurement

MEASUREMENT

Select the analysis **011 Cya.Ac. : 10 - 200mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT302.

I-020-p-Alkalinity TA: 2.0- 50.0° F

Reagent kit: 1MT134
Preparation time: ~ 6min



REAGENTS

Alkaphot P Tablets 1AP251

EQUIPMENT

Graduated Plastic Tube 14TP00
Crushing Rod 1AP018
Glass Tube 1CR099
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Alkaphot P tablet, crush it with the crushing rod and shake till it is dissolved.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 4 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **020 TA: 2.00 - 50.0 °F**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT045.

I-040- Aluminum: 0.05 - 3.00 mg/l

Reagent kit: 1MT136
Preparation time: ~ 5min

Photopod

LS

REAGENTS

AluminumbufferReagent	1RA010
Aluminum 1Reagent	1RA021
Alumimum 2Reagent	1RA030

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in a glass tube
Add 6 drops of Aluminumbuffer Reagent
Close andshake.
Add 6 drops of Aluminum 1Reagent
Close andshake.
Add 12 drops of Aluminum 2Reagent
Close andshake.
Wait 4 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **040 Al: 0.05 - 3.00 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT303.

I-041 - Aluminum: 0.20- 3.00 mg/L Al

Reagent kit: 1MT001
Preparation time: ~ 9min

Photopod

SP**REAGENTS**

Aluminum Tablets n°1 and n°2	1AP166
Demineralized water	1ED010

EQUIPMENT

Graduated Plastic Tube	14TP00*2
Glass Tube	1CR099
Crushing Rod	1AP018
Syringe 1 ml	1SU010
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 1 ml sample of water to analyze using the syringe, put it in the graduated plastic tube. Fill the Tube with demineralized water up to the 10 ml mark.

Put the cap and shake.

Add 1 Aluminum n°1 tablet, crush it with the crushing rod for 2 min and shake till it is dissolved. After complete dissolution of the tablet N°1, add 1 Aluminum n°2 tablet, crush it with the crushing rod and mix gently to dissolve by turning the tube. Mix by means of the rod to dissolve well the tablet and degas completely the mixture.

Fill a glass tube with this preparation using the plastic funnel then cover the tube.

Wait 5 minutes after crushing the Aluminum n°2 tablet (All the bubbles stemming from the effervescence must have disappeared).

Proceed to the measurement.

MEASUREMENT

Select the analysis **041 Al: 0.20- 3.00mg/L**

In another graduated plastic tube introduce 1 ml sample of water to analyze

Fill the Tube to the 10 ml mark with demineralized water.

Fill a glass tube with this preparation using the plastic funnel then cover the tube.

Put the black cover on top of the tube and press the key « zero »

Remove the tube and put the sample tube to analyze

Put the black cover on top of the tube and press the key « measure »

NOTA - Interference: Fluoruride and Polyphosphates

Please note that the refill kit's reference is 1MT304.

I-042- Aluminum: 0.02 - 0.30 mg/L Al

Reagent kit: 1MT001
Preparation time: ~ 8min

Photopod

SP**REAGENTS**

Aluminum Tablets n°1 and n°2 1AP166

EQUIPMENT

Graduated Plastic Tube 14TP00
Glass Tube 1CR099
Crushing Rod 1AP018
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Aluminum n°1 tablet, crush it with the crushing rod for 2 min and shake till it is dissolved.
After complete dissolution of the pill N°1, add 1 Aluminum n°2 tablet, crush it with the crushing rod and mix gently to dissolve by turning the tube. Mix by means of the rod to dissolve well the tablet and completely degas the mixture.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes after crushing the aluminum n°2 tablet (All the bubbles stemming from the effervescence must have disappeared).
Proceed to the measurement.

MEASUREMENT

Select the analysis **042 Al: 0.02 - 0.30 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA - Interference: Fluoruride and Polyphosphates

Please note that the refill kit's reference is 1MT304.

I-060–Ammonium: 0.08 - 1.60 mg/l NH₄-N
 METHOD COMPATIBLE WITH SEA WATER

Photopod
LS/SP

Reagent kit: 1MT193(or 1MT003 for sea water)
 Preparation time: ~ 11min

REAGENTS

Ammonia tablets 1 and 2	1AP152
Ammonia Conditioning Reagent	1AT170

EQUIPMENT

Graduated Plastic Tube	14TP00
Glass Tube	1CR099
Crushing Rod	1AP018
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.

If sample is sea water, add 1 spoonful of Ammonia Conditioning Reagent, shake to dissolve ~1min.

If turbidity appears, add 2 other spoonful of Ammonia Conditioning Reagent, shake to dissolve ~2 min.

Add the ammonia tablets 1 and 2, crush them with the crushing rod and shake to dissolve.

Fill a glass tube with this preparation using the plastic funnel then cover the tube.
 Wait 10 minutes.
 Proceed to measurement

MEASUREMENT

Select the analysis **060 NH₄-N: 0.08 -1.60mg/L**(Result in mg/L of N)
 Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
 Put the black cover on top of the tube and press the key « zero »
 Remove the tube and put the sample tube to analyze
 Put the black cover on top of the tube and press the key « measure »

To obtain the result as mg/L of NH₄⁺, multiply the result by 1,3

NB: If concentration is higher than 2 mg/l, the green color will be too strong and measurement won't be done. In this case, proceed to a dilution of the sample before adding reagents.

*Please note that the refill kit's reference is 1MT306.
 For sea water, the reference is 1MT358.*

I-061 - Ammonium: 0.20 - 4.80 mg/L N-NH₄

Reagent kit: 1MT002
Preparation time: ~ 6 min

Photopod

LS*REAGENTS*

Seignette Salt Reagent	1SD010
Nessler Reagent	1RD002

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the glass tube
Add 6 drops of Seignette Salt Reagent
Close and shake.
Add 6 drops of Nessler Reagent
Close and shake.
Wait 5 minutes
Proceed to the measurement

MEASUREMENT

Select the analysis **061 NH₄-N: 0.30 - 4.80 mg/L (Result in mg/L of N)**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

To obtain the result as mg/L of NH₄⁺, multiply the result by 1,3

Please note that the refill kit's reference is 1MT305.

I-062 –Ammonium: 0.80 - 24.0 mg/l N-NH₄

Reagent kit: 1MT002
Preparation time: ~ 6 min

Photopod

LS**REAGENTS**

Seignette Salt Reagent	1SD010
Nessler Reagent	1RD002

EQUIPMENT

Graduated Plastic Tube	14TP00
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 12,5 ml sample of water to analyze in the graduated plastic tube.
Add 8 drops of Seignette Salt
Close the tube and shake.
Add 8 drops of Nessler Reagent
Close the tube and shake.
Wait 5 minutes
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT


Select the analysis **062 NH₄-N: 1.0 - 24 mg/L** (*Result in mg/L of N*)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

To obtain the result as mg/L of NH₄⁺, multiply the result by 1,3

Please note that the refill kit's reference is 1MT305.

I-063-AmmoniumLR: 0,02-5 mg/L NH₄⁺-N

Reagent Kit: FTI535600
Preparation Time: ~ 6 min

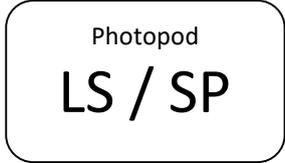

REAGENTS

Ammonia LR NH₄⁺ tubes
Vario Ammonia Salicylate F5 Powder Pack
Vario Ammonia Cyanurate F5 Powder pack
Demineralized water 1ED010

RECOMMENDED EQUIPMENT (consult us)

Automatic Pipette 1 - 5 ml 1PA023
Pipette Tip 1 - 5 ml 1EU013
24 tubes stand Ø16 1PT013

TEST INSTRUCTIONS

Take 2 Ammonia LR NH₄⁺tubes: one for the blank, the other for the sample.
With the pipette, put 2 ml of demineralized water in the blank tube.
Likewise, put 2 ml of water to analyse in the sample tube.
Add the content of oneVario Ammonia Salicylate F5 Powder pack in each tube then add the content of oneVario Ammonia Cyanurate F5 Powder packin each tube, close, and shake the tube for 30 seconds.
Wait 10 minutes.
Proceed to the measurement.

MEASUREMENT

Select the analysis **063 NH₄-N: 0,02 - 5 mg/L**
Take the tube for the blank and insert it in the photometer.
Put the black cover on top of the tube and press the key « zero ».
Remove the tube and put the sample tube to analyze.
Put the black cover on top of the tube and press the key « measure ».

To have the result as mg/l NH₄⁺, multiply the result by 1,3


INTERFERENCE

For strong alkaline or acidic water, you must adjust the pH at 7 by using hydrochloric acid 1 mol/L (if pH>7) or sodium hydroxide 1mol/L (if pH<7).
Iron interferes with the test.

I-064 - Ammonium HR: 0,5-50 mg/L NH₄⁺-N

Reagent Kit : FTI535650

Temps de préparation : ~ 6 min

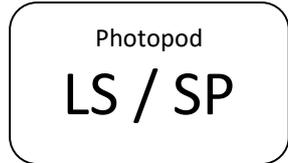**REAGENTS**Ammonia HR NH₄⁺ tubes

Vario Ammonia Salicylate F5 Powder Pack

Vario Ammonia Cyanurate F5 Powder pack

Demineralized water 1ED010

RECOMMENDED EQUIPMENT (consult us)

Automatic pipette 0,1 - 1 ml 1PA022

Pipette Tip 0,1 - 1 ml 1EU012

24 tubes stand Ø16 1PT013

TEST INSTRUCTIONSTake 2 Ammonia HR NH₄⁺tubes: one for the blank, the other for the sample.

With the pipette, put 0,1 ml of demineralized water in the blank tube.

Likewise, put 0,1 ml of water to analyse in the sample tube.

Add the content of oneVario Ammonia Salicylate F5 Powder pack in each tube then add the content of oneVario Ammonia Cyanurate F5 Powder pack in each tube, close, and shake the tube for 30 seconds.

Wait 10 minutes.

Proceed to the measurement.

MEASUREMENTSelect the analysis **064 NH₄-N: 0,5 - 50 mg/L**

Take the tube for the blank and insert it in the photometer.

Put the black cover on top of the tube and press the key « zero ».

Remove the tube and put the sample tube to analyze.

Put the black cover on top of the tube and press the key « measure ».

To have the result as mg/l NH₄⁺, multiply the result by 1,3**INTERFERENCE**

For strong alkaline or acidic water, you must adjust the pH at 7 by using hydrochloric acid 1 mol/L (if pH>7) or sodium hydroxide 1mol/L (if pH<7).

Iron interferes with the test.

I-072-Total Nitrogen LR: 0,3-20 mg/L N

Reagents kit reference: FTI2420703
Preparation time: ~ 80 min

Photopod
LS / SP

REAGENTS

Digestion tubes (empty tubes)
Blank tube (red label)
Reaction tubes
Digestion Reagent
Compensation Reagent
Nitrate-111

RECOMMENDED EQUIPMENT (consult us)

Automatic pipette 0,1 - 1 ml	1PA022
Pipette Tip 0,1 - 1 ml	1EU012
Automatic Pipette 1 - 5 ml	1PA023
Pipette Tip 1 - 5 ml	1EU013
24 tubes stand Ø16	1PT013
Wooden clamp	1PT007
Heating reactor	1RD010

TEST INSTRUCTIONS

Turn on the heating reactor. Preheat at 100 °C.
Take one Digestion tube.
With the pipette, put 5 ml of water to analyse in the tube.
Add 1 level scoop of No. 8 (black) digestion reagent, close, and shake for 30 seconds.
Put the tube in the reactor for 60 minutes at 100°C.
Remove the tube from the thermoreactor with the wooden clamp. (CAUTION: the tubes are hot!). Place the tube in the tube stand and allow to cool to room temperature.
Add 1 level scoop of No. 4 (white) compensation reagent, close and shake for 30 seconds.
This is the pre-treated sample.
Take 1 Reaction tube, and with the pipette, add 0,5 ml of pre-treated sample, close and return several times to mix the content (Caution: tube becomes warm!).
With the pipette, add 0,2 ml of Nitrate-111, close, and shake the tube.
Wait 10 minutes.

MEASUREMENT

Select the analysis **072 N: 0,3-20 mg/L**
Take the tube for the blank (tube with red label) and insert it in the photometer.
Put the black cover on top of the tube and press the key « zero ».
Remove the tube and put the sample tube to analyze.
Put the black cover on top of the tube and press the key « measure ».

NOTES

This test determines the inorganic compounds Ammonia, Nitrate and Nitrite, as well as organic compounds like amino acid, urea, complexing agents etc. Nitrogen compounds which are hardly to oxidise, as may be found in industrial sewage, are not digested or only partially.

I-073-Total Nitrogen HR: 3-200 mg/L N

Reagents kit reference: FTI2420703

Preparation time: ~ 80 min

Photopod

LS / SP**REAGENTS**

Digestion tubes (empty tubes)

Blank tube (red label)

Reaction tubes

Digestion Reagent

Compensation Reagent

Nitrate-111

Demineralized water

RECOMMENDED EQUIPMENT (consult us)

Automatic pipette 0,1 - 1 ml	1PA022
Pipette Tip 0,1 - 1 ml	1EU012
Automatic Pipette 1 - 5 ml	1PA023
Pipette Tip 1 - 5 ml	1EU013
24 tubes stand Ø16	1PT013
Wooden clamp	1PT007
Heating reactor	1RD010

TEST INSTRUCTION

Turn on the heating reactor. Preheat to 100 °C.

Take one Digestion tubes.

With the pipette, put 0,5 ml of water to analyse and, with the pipette 4.5 ml of demineralized water in the tube.

Add 1 level scoop of No. 8 (black) digestion reagent, close, and shake for 30 seconds.

Put the tube in the reactor for 60 minutes at 100°C.

Remove the tube from the thermoreactor with wooden clamp. (CAUTION: the tubes are hot!).

Place the tube in the tube stand and allow cooling to room temperature.

Add 1 level scoop of No. 4 (white) compensation reagent, close and shake for 30 seconds.

This is the pre-treated sample.

Take 1 reaction tube, and with the pipette, add 0,5 ml of pre-treated sample, close and return several times to mix the content (Caution: Tube becomes warm!).

With the pipette, add 0,2 ml of Nitrate-111, close, and mix the tube.

Wait 10 minutes.

MEASUREMENTSelect the analysis **073 N: 3-200 mg/L**

Take the tube for the blank (tube with red label) and insert it in the photometer.

Put the black cover on top of the tube and press the key « zero ».

Remove the tube and put the sample tube to analyze.

Put the black cover on top of the tube and press the key « measure ».

NOTES

This test determines the inorganic compounds Ammonia, Nitrate and Nitrite, as well as organic compounds like amino acid, urea, complexing agents etc. Nitrogen compounds

which are hardly to oxidise, as may be found in industrial sewage, are not digested or only partially.

I-080- Benzotriazole: 1.00 - 16.0 mg/L

Reagent kit: 1MT078
Preparation time: ~ 5.5min

Photopod

LS

REAGENTS

Triazole Reagent 1RT018

EQUIPMENT (consult us)

Graduated Glass Flask 125 ml	1FG000
UVLamp+ UV Protection Goggles	14LU01
UV Protection Goggles	FHA2113400
Glass Tube	1CR099
Plastic funnel	1EP021

TEST INSTRUCTIONS

Caution: the lamp produces UV rays hazardous to eyes and skin.
Wear UV protection goggles when the light is on.
Avoid touching the surface of the quartz of the lamp. Wipe the lamp after each use.

Take a 25 ml sample of water to analyze in the graduated glass flask.
Add 20 drops of Triazole Reagent (1RT018) and shake.
Put the UV glasses on.
Introduce the UV lamp in the flask and put it on for 5 minutes then turn it off.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement.

MEASUREMENT

Select the analysis **080 BZT: 1.00- 16.0 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA

To check that the lamp works correctly, take a solution at 5,0 mg/l of benzotriazole and make an analysis. If the result is below 5,0 mg/l, then change the lamp.

Please note that the refill kit's reference is 1MT307

I-101-Bromine: 1.00-13.5 mg/L Br₂

Reagent kit: 1MT138
Preparation time: ~ 7min



REAGENTS

DPD 1 Tablet	1D1018P
DPD Glycine Tablet	1NP000

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Total bromine test

Take a 12,5 ml sample of water to analyze in the graduated plastic tube.

In the presence of chlorine: Add 1 DPD Glycine tablet, crush it with the crushing rod and stir with the rod to dissolve.

Add 1 DPD 1 tablet, crush it with the crushing rod and stir with the rod to dissolve.

Fill a glass tube with this preparation using the plastic funnel then cover the tube.

Wait 5 minutes after crushing the tablet.

Proceed to the measurement

MEASUREMENT

Select the analysis **101 Br₂:1.00-13.5mg/L**

Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer

Put the black cover on top of the tube and press the key « zero »

Remove the tube and put the sample tube to analyze

Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT004.

I-102 - Bromine: 0.10 - 2.25 mg /L Br₂

Reagent kit: 1MT138
Preparation time: ~ 4min



REAGENTS

DPD 1 Tablet	1D1018P
DPD Glycine Tablet	1NP000

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Total brominetest

Take a 20 ml sample of water to analyze in the graduated plastic tube.
In the presence of chlorine: Add 1 DPD Glycine tablet; crush it with the crushing rod and stir with the rod to dissolve .
Add 1 DPD 1 tablet, crush it and stir with the rod to dissolve.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 2 minutes after crushing the tablet.
Proceed to the measurement

MEASUREMENT

Select the analysis **102 Br₂: 0.10 - 2.25 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT004.

I-110- Calcium: 20 - 200 mg/l CaCO₃

Reagent kit: 1MT139
Preparation time: ~ 4min

Photopod

SP

REAGENTS

Calcicol 1 Tablet	1AP252
Calcicol 2 Tablet	1AP252
Demineralized water	1ED010

EQUIPMENT

Graduated Plastic Tube	14TP00 x 2
Crushing Rod	1AP018
Glass Tube	1CR099
pH indicator 0 - 14	1PI110
Syringe 1 ml	1SU010
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take 1 ml sample of water to analyze using the syringe, put it in the graduated plastic tube. Fill the Tube to the 10 ml mark with Demineralized water. Put the cap and shake. Check with the pH indicator test strip that the pH is between 4 and 10 otherwise adjust it. Add 1 Calcicol 1 tablet, crush it with crushing rod and stir until dissolved ~ 30s. Add 1 Calcicol 2 tablet, crush it with crushing rod and stir until dissolved ~ 30s. Wait 2 minutes. Fill a glass tube with this preparation using the plastic funnel then cover the tube. Proceed to the measurement.

MEASUREMENT

Select the analysis **110 Ca: 20 - 200mg/L CaCO₃**
In a graduated plastic tube, introduce a 1 ml sample of water to analyze. Fill the tube to the 10 ml mark with Demineralized water. Fill a glass tube with this preparation using the plastic funnel then cover the tube. Put the black cover on top of the tube and press the key « zero ». Remove the tube and put the sample tube to analyze. Put the black cover on top of the tube and press the key « measure ».

NOTA: Concentration Ca mg/L = reading x 0.4

Interference: Mg < 200mg/L: nothing
Iron > 10 mg/L: concentration lower
Zinc > 5 mg/L: concentration higher

Please note that the refill kit's reference is 1MT309.

I-111 - Calcium: 2.0 -20.0mg/ICaCO₃

Reagent kit: 1MT139
Preparation time: ~ 3min



REAGENTS

Calcicol 1 Tablet	1AP252
Calcicol 2 Tablet	1AP252

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
pH indicator pH 0 - 14	1PI110
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Check with the pH indicator test strip that the pH is between 4 and 10 otherwise adjust it.
Add 1 Calcicol 1 tablet, crush it with crushing rod and stir until dissolved ~ 30s
Add 1 Calcicol 2 tablet, crush it with crushing rod and stir until dissolved ~ 30s
Wait 2 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **111 Ca.: 2.0- 20.0mg/L CaCO₃**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »


NOTA: Concentration Ca mg/L= reading x 0.4

Interference: Mg < 200 mg/L: nothing
Iron > 10 mg/L: concentration lower
Zinc > 5 mg/L: concentration higher

Please note that the refill kit's reference is 1MT309.

I-121 - Free chlorine and total chlorine: 0.50- 6.00 mg/L Cl₂

Reagent kit: 1MT140 and 1MT192
Preparation time: ~ 7min

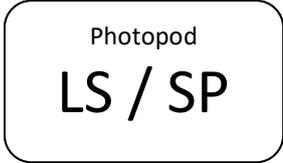

REAGENTS

DPD 1 Tablet (Free chlorine)	1D1018P
DPD 4 Tablet (total chlorine)	1D4004P

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (consult us)

TEST INSTRUCTIONS

Free chlorine test

Take a 12,5 ml sample of water to analyze in the graduated plastic tube.
Add 1 DPD 1 tablet, crush it and shake to dissolve
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes.
Proceed to the measurement.

Total chlorine test

Take a 12,5 ml sample of water to analyze in the graduated plastic tube.
Add 1 DPD 4 tablet, crush it and shake to dissolve
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes.
Proceed to the measurement.

MEASUREMENT

Select the analysis **121 Cl₂:0.50- 6.00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT116 (for free chlorine) & 1MT007 (for total chlorine).

I-122 -Free chlorine and total chlorine: 0.05- 1.00 mg /L

Reagent kit: 1MT140 and 1MT192

Preparation time: ~ 4min

Photopod

LS / SP

REAGENTS

DPD 1 Tablet (Free chlorine)	1D1018P
DPD 4 Tablet (total chlorine)	1D4004P

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Free chlorine test

Take a 20 ml sample of water to analyze in the graduated plastic tube.
Add 1 DPD 1 tablet, crush it and shake to dissolve
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 2 minutes.
Proceed to the measurement.

Total chlorine test

Take a 20 ml sample of water to analyze in the graduated plastic tube.
Add 1 DPD 4 tablet, crush it and shake to dissolve
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 2 minutes.
Proceed to the measurement.

MEASUREMENT

Select the analysis **122 Cl₂: 0.05 - 1.00 mg/L**

Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer

Put the black cover on top of the tube and press the key « zero »

Remove the tube and put the sample tube to analyze

Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT116 (for free chlorine) & 1MT007 (for total chlorine).

I-130-Chloride: 10 - 500 mg/L Cl⁻

Reagent kit: 1MT044
Preparation time: ~ 5min

Photopod

LS

REAGENTS

Chloride Reagent 1	1RC040
Chloride Reagent 2	1RC050
Demineralized water	1ED010

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013
Syringe 1 ml	1SU010

TEST INSTRUCTIONS

With the 10 ml syringe, introduce 9 ml of demineralized water in a glass tube
With a 1 ml syringe, introduce a 1 ml sample of water to analyze in the glass tube
Close the tube and shake.
Add 16 drops of Chloride Reagent 1
Close the tube and shake.
Add 16 drops of Chloride Reagent 2
Close the tube and shake.
Wait 3 minutes, invert the tube once every minute to homogenize.
Proceed to the measurement.

MEASUREMENT

Select the analysis **130 Cl⁻ : 10 - 500mg/L**
With the 10 ml syringe, introduce 9 ml of demineralized water in another glass tube
With the 1 ml syringe, introduce a 1 ml sample of water to analyze in the glass tube, put the cap, shake it and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA: Concentration in French degrees (°F) = reading x 0.14

Please note that the refill kit's reference is 1MT310.

I-131 - Chloride: 1.0 - 50.0 mg/L Cl⁻

Reagent kit: 1MT044
Preparation time: ~ 4min

Photopod

LS

REAGENTS

Chloride Reagent 1	1RC040
Chloride Reagent 2	1RC050

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the glass tube
Close glass tube and shake.
Add 16 drops of Chloride Reagent 1
Close glass tube and shake.
Add 16 drops of Chloride Reagent 2
Close glass tube and shake.
Wait 3 minutes, invert the tube once every minute to homogenize.
Proceed to the measurement.

MEASUREMENT

Select the analysis **131 Cl⁻ : 1.0 -50.0 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA: Concentration in French degrees (°F) = reading x 0.14

Please note that the refill kit's reference is 1MT310.

I-132 -Chloride: 5 - 200 mg/L Cl⁻

Reagent kit: 1MT141
Preparation time: ~ 5min

Photopod

SP

REAGENTS

Acidifying CD tablet	1AP268
Chloridol tablet	1AP268
Demineralized water	1ED010

EQUIPMENT

Graduated Plastic Tube	14TP00*2
Crushing Rod	1AP018
Glass Tube	1CR099
Syringe 1 ml	1SU010
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 1 ml sample of water to analyze using the 1ml syringe, put it in the graduated plastic tube. Fill the Tube up to the 10 ml mark with Demineralized water.

Close the tube and shake.

Add 1 acidifying CD tablet, crush it with the crushing rod and shake till it is dissolved ~ 1 min

Add 1 Chloridol tablet, wait 2 minutes, then crush it with the crushing rod and shake till it is dissolved ~ 30 s

Fill a glass tube with this preparation using the plastic funnel then cover the tube.

Proceed to the measurement

MEASUREMENT

Select the analysis **132 Cl⁻ : 5 - 200 mg/L**

In the graduated plastic tube introduce 1 ml sample of water to analyze

Fill the Tube to the 10 ml mark with Demineralized water.

Fill a glass tube with this preparation using the plastic funnel then cover the tube.

Put the black cover on top of the tube and press the key « zero »

Remove the tube and put the sample tube to analyze

Put the black cover on top of the tube and press the key « measure »

NOTA: Concentration in French degrees (°F) = reading x 0.14

Please note that the refill kit's reference is 1MT311.

I-133 -Chloride: 0.50-20.0 mg/L Cl⁻

Reagent kit: 1MT141
Preparation time: ~ 4min

Photopod

SP

REAGENTS

Acidifying CD tablet	1AP268
Chloridol tablet	1AP268
Demineralized water	1ED010

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 acidifying CD tablet, crush it with the crushing rod and shake till it is dissolved ~ 1 min
Add 1 Chloridol tablet, wait 2 minutes, then crush it with the crushing rod and shake till it is dissolved ~ 30 s
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **133 Cl⁻ :0.50-20.0 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA: Concentration in French degrees (°F) = reading x 0.14

Please note that the refill kit's reference is 1MT311.

I-140- Chromium VI: 0.10- 4.00 mg/L Cr⁶

Reagent kit: 1MT180
Preparation time: ~ 1.5min

Photopod

LS

REAGENTS

Chromium 1 Reagent	1RC032
Chromium 2 Reagent	1RC033
Sodium Fluoride Reagent	1SF000

EQUIPMENT

Glass Tube	1CR099
Syringe 5 ml	1SU012

TEST INSTRUCTIONS

With the syringe, take a 5 ml sample of water to analyze in the glass tube.
If water contains more than 1 mg/L of iron, eliminate it adding 3 drops of Sodium Fluoride
Add 4 drops of Chromium reagent 1 and shake.
Add 5 drops of Chromium reagent 2 and shake.
Wait 1 minute.
Proceed to the measurement.

MEASUREMENT

Select the analysis **140 Cr6:0.10 - 4.00 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA: Chromium VI in mg/l CrO₄²⁻ = reading x 2,23

Please note that the refill kit's reference is 1MT009.

I-141 - Chromium VI: 0.05 - 2.00 mg/L Cr⁶

Reagent kit: 1MT142
Preparation time: ~ 6min

Photopod

SP

REAGENTS

Chromicol 1 tablet	1AP281
Chromicol 2 tablet	1AP281

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Chromicol 1 tablet, crush it with the crushing rod and shake till it is dissolved ~ 30s
Add 1 Chromicol 2 tablet, crush it with the crushing rod and shake till it is dissolved ~ 30s
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **141 Cr6:0,05 - 2,00 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

INTERFERENCES

Dissolved iron concentration greater than 1 mg / L will give lower results in chrome.
To increase sensitivity then add 2 chromicol 1 tablet and 1 chromicol 2 tablet.
This instruction may not be applied if the sample matrix contains tannin.

Please note that the refill kit's reference is 1MT312.

I-150-FreeCopper: 0.05 - 5.0 mg/L Cu

Reagent kit: 1MT181
Preparation time: ~ 3.5min

Photopod

LS

REAGENTS

Copper 1 Reagent	1RC060
Copper 2 Reagent	1RC070

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

With the syringe, take a 10 ml sample of water to analyze and introduce in a glass tube
Add 5 drops of Copper 1 Reagent
Close and shake
Add 5 drops of Copper 2 Reagent
Close and shake
Wait 3 minutes
Proceed to the measurement

MEASUREMENT

Select the analysis **150 Cu: 0.05 - 5.00 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT313.

**I-151-Free Copper, Total Copper and Copper chelated:
0,20 - 5,00 mg/L Cu**

Photopod

SP

Reagent kit: 1MT011
Preparation time: ~ 6min

REAGENTS

Copper n°1 Tablet	1AP186
Copper n°2 Tablet	1AP186

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Free Copper

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Copper n°1 tablet. Crush it with the crushing rod and shake to dissolve. ~30 s
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes after crushing the tablet.
Proceed to the measurement.

Total Copper

Transfer the content of the glass tube in the graduated plastic tube.
Add 1 Copper n°2 tablet. Crush it with the crushing rod and shake to dissolve. ~15 s
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement.

Chelated Copper

Chelated Copper = Total Copper - Free Copper

MEASUREMENT

Select the analysis **151 Cu: 0.20- 5.00 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT314.

I-160- Cyanide: 0,02 - 0,50 mg/L CN⁻

Reagent kit: 1MT012
Preparation time: ~ 11min

Photopod

LS

REAGENTS

Cyanide n°1 Reagent	1RC005
Cyanide n°2 Reagent	1RC008
Cyanide n°3 Reagent	1RC011
Cyanide buffer Reagent	1RC018

EQUIPMENT

Graduated Plastic Tube	14TP00
Plastic spoon	1J0000
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyse in the graduated plastic tube.
Add 4 drops of Cyanide buffer Reagent and shake
Add 4 drops of Cyanide n°1 Reagent and shake
Wait 1 minute.
Add 1 plastic spoon to the brim of Cyanide n°2 Reagent and shake
Wait 2 minutes.
Add 16 drops of Cyanide n°3 Reagent and shake
Wait 7 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement.


MEASUREMENT

Select the analysis **160 CN- : 0,02 - 0,50 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT315.

I-174 - COD HR: 0,5 - 15 g/L O₂ (500 - 15000 mg/L O₂)

Reagents kit reference: 12DC02
Preparation time: ~ 2h30

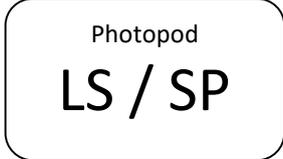

ATTENTION: The tubes contains sulfuric acid <90% (corrosive) and potassium dichromate <0,5% (toxic). Before starting measure, please read MSDS.

REAGENTS

Reaction tubes
Demineralized water 1ED010

RECOMMENDED EQUIPMENT

Automatic pipette 0,1 - 1 ml 1PA022
Pipette Tip 0,1 - 1 ml 1EU012
24 tubes stand Ø16 1PT013
Wooden clamp 1PT007
Heating reactor 1RD010

TEST INSTRUCTIONS (consult us)

Turn on the heating reactor. Preheat to 150 °C.
Take 2 tubes: one for the blank and one for the sample.
With the pipette, put 0,2 ml of demineralized water in the tube for the blank, close and shake gently the tube. Likewise, put 0,2 ml of water to analyze in the tube for the sample.

Be careful, the reaction is exothermic and the tube becomes hot.

Put the tubes in the heating reactor at 150°C during 120 minutes.
After 120 minutes, take the tubes with the wooden clamp (be careful, they are very hot) and shake gently. Put the tube in the tube stand et let cool to room temperature (>20 minutes).

MEASUREMENT

Select the analysis **174DCO: 0,5-15g/L**
Take the tube for the blank tube and insert it in the photometer.
Put the black cover on top of the tube and press the key « zero ».
Remove the tube and put the sample tube to analyze.
Put the black cover on top of the tube and press the key « measure ».


INTERFERENCES

Suspended solids in the tube lead to incorrect measurements. For this reason it is important to place the tubes carefully in the sample chamber. The precipitate at the bottom of the sample should be not suspended. Samples can be measured when the Chloride content does not exceed 1000 mg/l.
If the sample is COD < 1g/L, we advise to do measure with kit COD MR for a best precision.
If the sample is COD < 0,1g/L, we advise to do measure with kit COD LR for a best precision.

I-175-COD MR: 50 - 1500 mg/L O₂

Reagents kit reference: 12DC01

Preparation time: ~ 2h30

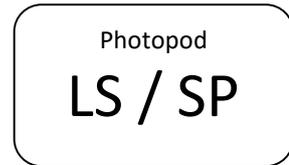

ATTENTION: The tubes contains sulfuric acid <90% (corrosive) and potassium dichromate <0,5% (toxic). Before starting measure, please read MSDS.

REAGENTS

Reaction tubes

Demineralized water

RECOMMENDED EQUIPMENT (consult us)

Automatic Pipette 1 - 5 ml	1PA023
Pipette Tip 1 - 5 ml	1EU013
24 tubes stand Ø16	1PT013
Wooden clamp	1PT007
Heating reactor	1RD010

TEST INSTRUCTIONS

Turn on the heating reactor. Preheat to 150 °C.

Take 2 tubes: one for the blank and one for the sample.

With the pipette, put 2 ml of demineralized water in the tube for the blank, close and shake gently the tube. Likewise, put 2 ml of water to analyze in the tube for the sample.

Be careful, the reaction is exothermic and the tube becomes hot.

Put the tubes in the heating reactor at 150°C during 120 minutes.

After 120 minutes, take the tubes with the wooden clamp (be careful, they are very hot) and shake gently. Put the tubes in the tube stand et let cool to room temperature (>20 minutes).

MEASUREMENT

Select the analysis **175 DCO: 50-1500 mg/L**

Take the tube for the blank tube and insert it in the photometer.

Put the black cover on top of the tube and press the key « zero ».

Remove the tube and put the sample tube to analyze.

Put the black cover on top of the tube and press the key « measure ».

INTERFERENCES

Suspended solids in the tube lead to incorrect measurements. For this reason it is important to place the tubes carefully in the sample chamber. The precipitate at the bottom of the sample should be not suspended. Samples can be measured when the Chloride content does not exceed 1000 mg/l.

If the sample is COD < 100 mg/L, we advise to do measure with kit COD LR for a best precision.

I-176-COD LR: 5 - 150 mg/L O₂

Reagents kit reference: 12DC00 ou FTI2420720
Preparation time: ~ 2h30

Photopod
LS / SP

ATTENTION: The tubes contains sulfuric acid <90% (corrosive) and potassium dichromate < 0,5% (toxic). Before starting measure, please read MSDS.

REAGENTS

Reaction tubes
Demineralized water 1ED010

RECOMMENDED EQUIPMENT

Automatic Pipette 1 - 5 ml 1PA023
Pipette Tip 1 - 5 ml 1EU013
24 tubes stand Ø16 1PT013
Wooden clamp 1PT007
Heating reactor 1RD010

TEST INSTRUCTION

Turn on the heating reactor. Preheat to 150 °C.
Take 2 tubes: one for the blank and one for the sample.
With the pipette, put 2 ml of demineralized water in the tube for the blank, close and shake gently the tube. Likewise, put 2 ml of water to analyze in the tube for the sample.

Be careful, the reaction is exothermic and the tube becomes hot.

Put the tubes in the heating reactor at 150°C during 120 minutes.
After 120 minutes, remove the tubes with the wooden clamp (be careful, they are very hot) and shake gently. Put the tubes in the tube stand et allow to cool to room temperature (>20 minutes).

MEASUREMENT

Select the analysis **176 DCO: 5-150 mg/L**
Beware; if Photopod is used with Odeon for this method, sample tube has to be put in the device before blank tube.
Insert **the sample tube** in the photometer, put the black cover and press the key « zero ».
Remove the tube and put **the blank tube**
Put the black cover on top of the tube and press the key « measure ».

If Photopod is used with Spectralab, put the blank tube first, then the sample tube.

INTERFERENCE

Suspended solids in the tube lead to incorrect measurements. For this reason it is important to place the tubes carefully in the sample chamber. The precipitate at the bottom of the sample should be not suspended. Samples can be measured when the Chloride content does not exceed 1000 mg/l.

I-180- DEHA : 0.02 - 1.00 mg/L

Reagent kit: 1MT182
Preparation time: ~ 11.5min

Photopod

LS

REAGENTS

DEHA n°1 Reagent	1RD012
DEHA n°2 Reagent	1RD013

EQUIPMENT

Graduated Plastic Tube	14TP00
Glass Tube	1CR099
Spoon	1J0000
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 spoonful of DEHA n°1 reagent and shake strongly ~30 s
Add 5 drops of DEHA n°2 Reagent, and shake
Wait 10 minutes, inverting the tube every minute to homogenize.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **180 DEHA: 0,02 - 1,00 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA- Avoid exposure to sunlight.
Make the measurement at a temperature between 22°C and 28°C.

INTERFERENCES


Reagents react with iron. Presence of Iron will give overestimated results

substance	concentration	substance	concentration
Borate	500 mg/l	Nickel	0,8 mg/l
Cobalt	0,025 mg/l	Phosphate	10 mg/l
Copper	8,0 mg/l	Phosphonates	10 mg/l
Hardness	1000 mg/l	Sulfate	1000 mg/l
Lignosulfonates	0,05 mg/l	Zinc	50 mg/l
Manganese	0,8 mg/l		
Molybdene	80 mg/l		

Please note that the refill kit's reference is 1MT112.

I-191 - Chlorine Dioxide: 2.4-28.5mg/L ClO₂

Reagent kit: 1MT177
Preparation time: ~ 7min

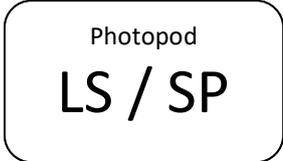

REAGENTS

DPD 1 Tablet	1D1018P
Glycine Tablet	1NP000

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 12.5 ml sample of water to analyze in the graduated plastic tube.
Add 1 Glycine tablet and shake to obtain total dissolution ~ 1min
Add 1 DPD 1 tablet and shake to obtain total dissolution ~ 1min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement 5 min after crushing the DPD 1 tablet.


MEASUREMENT

Select the analysis **191 ClO₂: 2.40-28.5 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT069.

I-192 - Chlorine Dioxide: 0.20 - 4.75 mg/L ClO₂

Reagent kit: 1MT177
Preparation time: ~ 4min

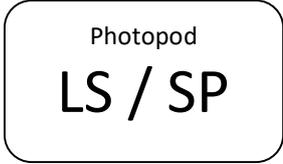

REAGENTS

DPD 1 Tablet	1D1018P
Glycine Tablet	1NP000

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 20 ml sample of water to analyze in the graduated plastic tube.
Add 1 Glycine tablet and shake to obtain total dissolution ~ 1 min
Add 1 DPD 1 tablet and shake to obtain total dissolution ~ 1 min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement 2 min after crushing the DPD 1 tablet.

MEASUREMENT

Select the analysis **192 ClO₂: 0.20 -4.75 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT069.

I-200-Total Hardness: 5.0 -50.0°F

Reagent kit: 1MT143
Preparation time: ~ 5min



REAGENTS

Hardicol n°1Tablet	1AP254
Hardicoln°2Tablet	1AP254
Demineralized water	1ED010

EQUIPMENT

Graduated Plastic Tube	14TP00* 2
Crushing Rod	1AP018
Glass Tube	1CR099
Syringe 5 ml	1SU012
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

With the syringe, take 4 ml of water to analyze, introduce it in the graduated plastic tube and complete with demineralized water up to 10 ml
Close the tube and shake.
Add 1 Hardicol n°1 tablet, crush it with the crushing rod and shake to dissolve ~ 1min
Add 1 Hardicol n°2 tablet, crush it with the crushing rod and shake to dissolve ~ 30s
Ensure that the tablet are well dissolved
Wait 2 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT


Select the analysis **200 TH: 5.0 -50.0 °F**
With the syringe, take 4 ml of water to analyze, introduce it in the graduated plastic tube and complete with demineralized water up to 10 ml
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA - For water containing Iron in concentrations higher than 10 mg/l, results will be underestimated.
The pH of the water should be between 4 and 10.

Please note that the refill kit's reference is 1MT047.

I-201- Total Hardness: 2.0 -20.0°F

Reagent kit: 1MT143
Preparation time: ~ 4min

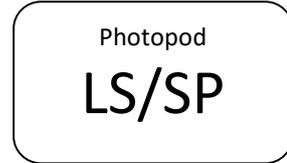

REAGENTS

Hardicol n°1Tablet	1AP254
Hardicol n°2Tablet	1AP254

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Hardicol n°1 tablet, crush it with the crushing rod and shake to dissolve~ 1 min
Add 1 Hardicol n°2 tablet, crush it with the crushing rod and shake to dissolve~ 30 s
Ensure tha the tablet are well dissolved
Wait 2 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **201 TH: 2.0 - 20.0 °F**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key« zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA - For water containing Iron in concentrations higher than 10 mg/l, results will be underestimated.
The pH of the water should be between 4 and 10.

Please note that the refill kit's reference is 1MT047.

I-210- Iron: 0,05 - 5.00 mg/L Fe²⁺ Fe³⁺

Reagent kit: 1MT144
Preparation time: ~ 3min

Photopod

LS

REAGENTS

FerrordisReagent	14F600
Chlorhydric Acid ½ (optional)	1AC000
Sodium Hydroxide 1N (optional)	1SH055

EQUIPMENT

pH indicator test strips 0-14	1PI110
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Check with the pH indicator test strip that the pH of the water to analyze is between 3 and 9 (T° ideally between 15 and 25°C), if not, adjust with Chlorhydric Acid or Sodium Hydroxide.
Take a 10 ml sample of water to analyze in a glass tube
Add 6 drops of Ferrordis reagent, shake.
Wait 2 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement.

MEASUREMENT

Select the analysis **210 Fe: 0,05- 5.00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT317.

I-211 - Iron: 0.2 - 20,0 mg/L Fe²⁺ Fe³⁺

Reagent kit: 1MT145
Preparation time: ~ 4min

Photopod

SP

REAGENTS

Iron 1 HR Tablet	1AP156
Chlorhydric Acid ½ (optional)	1AC000
Sodium Hydroxide 1N (optional)	1SH055

EQUIPMENT

pH indicator test strips 0-14	1PI110
Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Check with the pH indicator test strip that the pH of the water to analyze is between 3 and 9 (T° ideally between 15 and 25°C), if not, adjust with Chlorhydric Acid or Sodium Hydroxide.
Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Iron 1 HR tablet; crush it with the crushing rod.
Close the tube and shake strongly to dissolve ~ 3 min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 1 min
Proceed to the measurement.

MEASUREMENT

Select the analysis **211 Fe: 0,2-20.0 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT318.

I-212 - Iron: 0.05 - 5.00 mg/L Fe²⁺ Fe³⁺

Reagent kit: 1MT146
Preparation time: ~ 7min

Photopod

SP

REAGENTS

Iron MR 1 Tablet	1AP292
Iron MR 2 Tablet	1AP292

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Iron MR1 tablet, crush it with the crushing rod and shake to dissolve ~ 30 s
Add 1 Iron MR2 tablet, crush it with the crushing rod and shake to dissolve ~ 1 min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 min
Proceed to the measurement.

MEASUREMENT

Select the analysis **212 Fe: 0.05 -5.00 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT319

I-213- Iron: 0.1-10.0 mg/L Fe²⁺ Fe³⁺

Reagent kit: 1MT194
Preparation time: ~ 13 min

Photopod

LS

REAGENTS

Iron Reagent 1	1RF005
Iron Reagent 2	1RF006
Iron Reagent 3	1RF007

EQUIPMENT

Graduated Plastic Tube	14TP00
Glass Tube	1CR099
Plastic spoon	1J0000
Plastic funnel	1EP021 (<i>consult us</i>)
pH indicator paper 0-14	1PI110 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 20 ml sample of water to analyze in the graduated plastic tube.
Check with the pH indicator test strip that the pH of the water to analyze is between 3 and 9.
Add 10 drops of Iron reagent 1 and shake.
Add 1 spoonful of Iron reagent 2 and shake to dissolve.
Add 10 drops of Iron reagent 3 and shake.
Wait 10 minutes.
Fill a glass tube with this preparation by using the plastic funnel, then cover the tube.
Proceed to the measurement.

MEASUREMENT

Select the analysis **213 Fe: 0,1- 10,0 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA - To measure Iron II only, proceed the same way as Total Iron measurement but without adding Iron Reagent 2.

Please note that the refill kit's reference is 1MT359

I-220-Fluoride: 0.10- 2.00 mg/L F⁻

Reagent kit: 1MT110
Preparation time: ~ 5.5min



REAGENTS

Fluoride test tubes 14RF06

EQUIPMENT

Syringe 2 ml 1SU001
Glass Tube 1CR099

Option for more accurate analysis

Automatic Pipette 1 - 5 ml 1PA023
+Pipette Tips 1 - 5 ml 1EU003
Or
Graduated Pipette 2 ml 1PG001
+Macropipette 1T0007

TEST INSTRUCTIONS

Take a 2 ml sample of water to analyze and insert it in a Fluoride test tube, cover the tube and invert it 3 times to homogenize.
Wait 5 minutes.
Proceed to the measurement.

MEASUREMENT

Select the analysis **220 F⁻:0,10- 2.00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer.
Put the black cover on top of the tube and press the key « zero ».
Remove the tube and put the sample tube to analyze.
Put the black cover on top of the tube and press the key « measure ».

NOTA -Aluminum, calcium and iron disturb the reaction and lead to underestimated results.
Nitrates interfere when superior to 100 mg/l.

I-221 - Fluoride: 0.20- 2.00 mg/L F⁻

Reagent kit: 1MT147
Preparation time: ~ 7min

Photopod

SP

REAGENTS

Fluoride 1 Tablet	1AP179
Fluoride 2 Tablet	1AP179

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Fluoride 1 tablet, crush it with the crushing rod and stir with the rod to dissolve ~ 30 s
Don't shake the plastic tube
Add 1 Fluoride 2 tablet, crush it with the crushing rod and stir with the rod to dissolve ~ 90 s
Don't shake the plastic tube
Wait 5 min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement.

MEASUREMENT

Select the analysis **221 F⁻: 0,20- 2.00 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT320.

I-230- Hydrazine : 0,10-1,00 mg/L N₂H₄

Reagent kit: 1MT019
Preparation time: ~ 3min

Photopod

LS

REAGENTS

DAB Indicator 1ID001

EQUIPMENT

Graduated Plastic Tube 14TP00
Glass Tube 1CR099
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 5 ml sample of water to analyze in the graduated plastic tube.
Add DAB indicator to 10 ml and shake
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 2 minutes.
Proceed to the measurement


MEASUREMENT

Select the analysis **230 N₂H₄:0,10 - 1,00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT323.

I-240 -Magnesium: 5.0 - 50.0 mg/L Mg

Reagent kit: 1MT161
Preparation time: ~ 5min

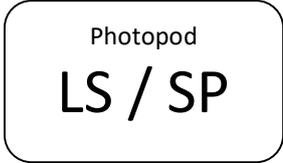

REAGENTS

MagnecolTablet	1AP193
Demineralized water	1ED010

EQUIPMENT

Graduated Plastic Tube	14TP00*2
Crushing Rod	1AP018
Glass Tube	1CR099
Syringe 1 ml	1SU010
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

With the syringe, take a 1 ml sample of water to analyze, introduce it in the graduated tube then complete up to 10 ml with demineralized water.
Close and shake
Add 1 Magnecol tablet; crush it with the crushing rod and shake to dissolve ~ 1min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 3 min
Proceed to the measurement


MEASUREMENT

Select the analysis **240 Mg: 5.00 - 50.0 mg/L**
In the graduated plastic tube introduce 1 ml sample of water to analyze
Fill the Tube to the 10 ml mark with Demineralized water.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT325.

I-241 - Magnesium: 0.50- 5.00 mg/L Mg

Reagent kit: 1MT161
Preparation time: ~ 4 min

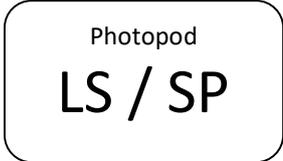

REAGENTS

Magnecol Tablet 1AP193

EQUIPMENT

Graduated Plastic Tube 14TP00
Crushing Rod 1AP018
Glass Tube 1CR099
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Magnecol tablet, crush it with the crushing rod and shake to dissolve ~ 1 min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 3 min
Proceed to the measurement

MEASUREMENT

Select the analysis **241 Mg:0.05-5.00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure ».

Please note that the refill kit's reference is 1MT325.

I-250 - Manganese: 0.20- 5.00 mg/L Mn

Reagent kit: 1MT050
Preparation time: ~ 6min

Photopod

LS

REAGENTS

Manganese 1 Reagent	1RM007
Manganese 2 Reagent	1RM008
Manganese 3 Reagent	1RM009

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the glass tube
Add 8 drops of Manganese 1 Reagent
Close and shake.
Add 8 drops of Manganese 2 Reagent
Close and shake.
Wait 2 minutes.
Add 8 drops of Manganese 3 Reagent
Close and shake.
Wait 5 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **250 Mn: 0,20 - 5,00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

INTERFERENCES

Concentrations in Ca^{2+} and Mg^{2+} higher than 300 mg/l lead to over-estimated results.
In the presence of Ca^{2+} , concentration in phosphates higher than 5 mg/l lead to underestimated results.
pH of the sample should be between 3 and 10. Temperature of the sample should be between 15 and 25°C.

Please note that the refill kit's reference is 1MT326.

I-251 - Manganese: 0.10- 8.00 mg/L Mn

Reagent kit: 1MT162
Preparation time: ~ 6min

Photopod

SP

REAGENTS

Manganese HR 1 Tablet	1AP174
Manganese HR 2 Tablet	1AP174

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Manganese HR 1 tablet, crush it with the crushing rod and shake to dissolve ~ 45s
Add 1 Manganese HR 2 tablet, crush it with the crushing rod and shake to dissolve ~ 45 s
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes
Proceed to the measurement

MEASUREMENT

Select the analysis **251 Mn: 0.10- 8.00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure ».

Please note that the refill kit's reference is 1MT327.

I-270 - Molybdates: 0.5 - 20.0 mg/L MoO₄-Mo

Reagent kit: 1MT183
Preparation time: ~ 1min

Photopod

LS

REAGENTS

Molybdate Reagent compensator 1RM010
Molybdate reagent 1RM016

EQUIPMENT

Glass Tube 1CR099
Syringe 10 ml 1SU007
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

With syringe, take a 10 ml sample of water to analyze in the glass tube
Add 5 drops of Molybdate Reagent compensator
Close and shake
Add 5 drops of Molybdate Reagent.
Close and shake
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 1 minute.
Proceed to the measurement

MEASUREMENT

Select the analysis **270 MoO₄-Mo: 0.5 - 20.0 mg/L** (*Result in mg/l de Mo*)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Concentration as MoO₄ mg/l = result x 1,66

Concentration as Na₂MoO₄ mg/l = result x 2,15

Please note that the refill kit's reference is 1MT329.

I-271 - Molybdates: 3.0 - 60.0 mg/L MoO₄-Mo

Reagent kit: 1MT024
Preparation time: ~ 2min

Photopod

SP

REAGENTS

Molybdate n°1 Tablet	1AP175
Molybdate n°2 Tablet	1AP175

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 molybdate n°1 tablet, crush it with the crushing rod and shake to dissolve ~ 1min
Add 1 molybdate n°2 tablet, crush it with the crushing rod and shake to dissolve ~ 1min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **271 MoO₄-Mo: 3.0- 60.0 mg/L** (*Result in mg/l of Mo*)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Concentration as MoO₄ mg/l = result x 1,66

Concentration as Na₂MoO₄ mg/l = result x 2,15

Please note that the refill kit's reference is 1MT330.

I-272 - Molybdates: 20 - 200 mg/L MoO₄-Mo

Reagent kit: 1MT183
Preparation time: ~ 1.5 min

Photopod

LS

REAGENTS

Molybdate Reagent compensator 1RM010
Molybdate reagent 1RM016

EQUIPMENT

Glass Tube 1CR099
Syringe 10 ml 1SU007
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

With the syringe, take a 10 ml sample of water to analyze in the glass tube
Add 5 drops of Molybdate Reagent compensator
Close and shake
Add 5 drops of Molybdate Reagent.
Close and shake
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 1 minute.
Proceed to the measurement

MEASUREMENT

Select the analysis **272 MoO₄-Mo: 20 - 200 mg/L** (*Result in mg/l of Mo*)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »


Concentration as MoO₄ mg/l = result x 1,66

Concentration as Na₂MoO₄ mg/l = result x 2,15

Please note that the refill kit's reference is 1MT329.

I-280- Nickel: 0,10-5,00 mg/L Ni

Reagent kit: 1MT164
Preparation time: ~ 4min

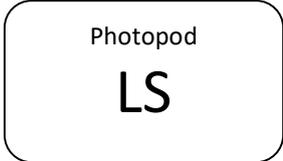

REAGENTS

Nickel 1 Reagent	1RN011
Nickel 2 Reagent	1RN012

EQUIPMENT

Graduated Plastic Tube	14TP00
Plastic Spoon	1J0000
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 spoonfull to the brim of of Nickel 1 Reagent and shake
Add 10 drops of Nickel 2 Reagent, and shake (color changes to orange).
Wait 3 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **280 Ni: 0,10 - 5,0mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA - Interferences happen when:

Mn ²⁺	>	1mg/l
Co ²⁺ Cu ²⁺ Fe ³⁺	>	5mg/l
Cr ³⁺ Zn ²⁺	>	10 mg/l

Please note that the refill kit's reference is 1MT331.

I-281 -Nickel: 0.50- 10 mg/L Ni

Reagent kit: 1MT079
Preparation time: ~3 min

Photopod

SP

REAGENTS

Nickeltest 1 tablet	1AP284
Nickeltest PR POWDER	1AP284
Nickeltest 2 tablet	1AP284

EQUIPMENT

Graduated Plastic Tube	14TP00
Plastic spoon	1J0000
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Nickeltest 1 tablet, crush it with the crushing rod and shake to dissolve. ~30 s
Only if the sample contains iron: Add 1 spoonful of Nickeltest PR POWDER and shake
Add 1 Nickeltest 2 tablet, crush it with the crushing rod and shake to dissolve ~30 s
A bit of foam appears. Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 2 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **281 Ni: 0.50- 10mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

NOTA - Interferences happen when:

Mn ²⁺	>	1mg/l
Co ²⁺ Cu ²⁺ Fe ³⁺	>	5mg/l
Cr ³⁺ Zn ²⁺	>	10 mg/l

Please note that the refill kit's reference is 1MT332.

I-300 - Nitrates : 0.10 - 1.00 mg/L NO₃⁻-N

Photopod

SP

Reagent kit: 1MT101
Preparation time: ~ 17min

REAGENTS

Nitrate Reagents	1AP163
Nitratetest Powder	1AP163
Nitratetest Tablets (contained in a white flask)	1AP163
Nitricol Tablets	1AP163

EQUIPMENT

Graduated Plastic Tube	14TP00 x 2
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 20 ml sample of water to analyze in the graduated plastic tube.
Add 1 spoonful of Nitratetest Powder and 1 Nitratetest tablet (from the white flask). Don't crush the tablet. Cover the tube and shake for 1 minute.
Wait 1 minute, then invert the tube 4 times to allow the flocculation. Wait until the liquid is clear (~2 minutes).
Remove the cap and wipe around the top of the tube with clean tissue.
Pour carefully 10 ml of the clear solution in another graduated plastic tube.
Add 1 Nitricol tablet, crush it with the crushing rod and shake to dissolve.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 10 minutes (after crushing the Nitricol tablet)
Proceed to the measurement.

MEASUREMENT

Select the analysis **300 NO₃-N: 0.10 -1.00 mg/L** (*Result in mg/L of N*)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer.
Put the black cover on top of the tube and press the key « zero ».
Remove the tube and put the sample tube to analyze.
Put the black cover on top of the tube and press the key « measure ».

Concentration as mg/L NO₃⁻ = result x 4,4

Nitrite Correction

The reagent also reacts with Nitrite. Most of the time, concentration in Nitrite is low compared to concentration in nitrate. But concentration in nitrite can be measured (in mg/l) and deduced from the concentration measured with this method.

Please note that the refill kit's reference is 1MT333.

I-301 - Nitrates : 0,06 - 2,30 mg/L NO₃⁻-N
 METHOD COMPATIBLE WITH SEA WATER



Reagent kit: 1MT184
 Preparation time: ~ 10 min

REAGENTS

Nitrates 1 Reagent	1RN015
Powder nitrate Reagent	1PN010

EQUIPMENT

Graduated Plastic Tube	14TP00 x 2
Crushing Rod	1AP018
Glass Tube	1CR099
Syringe 10ml	1SU013
Filter holder	14PF09
Filter paper	14PF05
Clamp	1PM010 (<i>consult us</i>)
Plastic spoon	1J0000

TEST INSTRUCTIONS

Take the filter holder, unscrew it and introduce it with the clamp, screw the bracket.
 In the first graduated plastic tube, take a 10 ml sample of water to analyze
 Add 1 spoonfull of Powder nitrate Reagent
 Close and shake 1 minute
 With the syringe 10ml take totally sample (water+zinc)
 Fit the 10 ml syringe containing the sample on the filter holder and gently squeeze out a few drops of (=rinsing the filter).
 Filter the sample and introduce it in the second plastic tube graduated to 5 ml graduation.
 Add nitrate 1 Reagent to 10 ml and shake.
 Wait 3 minutes.
 Proceed to the measurement
 Unscrew the filter holder, remove the soiled filter and clean all the water.

MEASUREMENT

Select the analysis **301 NO₃-N: 0,06 - 2.30 mg/L** (*Result in mg/L of N*)
 Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
 Put the black cover on top of the tube and press the key « zero »
 Remove the tube and put the sample tube to analyze
 Put the black cover on top of the tube and press the key « measure »

Concentration as mg/L NO₃⁻ = result x 4,4

Nitrite correction

The method optionally reacts with the nitrite present in the sample. Nitrite content is low compared to that of nitrate in most of the water, but can determine the concentration (mg / l N) of nitrite and then deduct the value of the content (in mg / l N) nitrate.

Please note that the refill kit's reference is 1MT350.

I-302 - Nitrates: 1.0 - 22.5 mg/L NO₃⁻-N
I-304 - Nitrates: 4,5 - 45.0 mg/L NO₃⁻-N

Reagent kit: 1MT101
 Preparation time: ~ 17 min



REAGENTS

Nitrate Reagents	1AP163
Nitratest Powder	1AP163
Nitratest Tablets (contained in a white flask)	1AP163
Nitricol Tablets	1AP163
Demineralized water	1ED010

EQUIPMENT

Graduated Plastic Tube	14TP00 x 3
Crushing Rod	1AP018
Syringe 1 ml	1SU010
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

With the syringe take a 1 ml sample of water to analyze, introduce it in the graduated plastic tube then complete up to 20 ml with Demineralized water
 Close and shake
 Add 1 spoonful of Nitratest Powder and 1 Nitratest tablet (from the white flask). Don't crush the tablet. Cover the tube and shake for 1 minute.
 Wait 1 minute, and then invert the tube 4 times to allow the flocculation. Wait until the liquid is clear (~2 minutes).
 Remove the cap and wipe around the top of the tube with a clean tissue.
 Pour carefully 10 ml of the clear solution in another graduated plastic tube.
 Add 1 Nitricol tablet, crush it with the crushing rod and shake to dissolve.
 Fill a glass tube with this preparation using the plastic funnel then cover the tube.
 Wait 10 minutes (after crushing the Nitricol tablet)
 Proceed to the measurement.

MEASUREMENT

Select the analysis **302 NO₃: 1,0 - 22,5 mg/L** (*Result in mg/L of N*)
 304 NO₃-N: 4,5 - 45.0 mg/L (*Result in mg/L of N*)

In the graduated plastic tube introduce 1 ml sample of water to analyze
 Fill the Tube to the 20 ml mark with Demineralized water.
 Fill a glass tube with this preparation using the plastic funnel then cover the tube.
 Put the black cover on top of the tube and press the key « zero »
 Remove the tube and put the sample tube to analyze
 Put the black cover on top of the tube and press the key « measure »

Concentration as mg/L NO₃⁻ = result x 4,4

Nitrite correction

The method optionally reacts with the nitrite present in the sample. Nitrite content is low compared to that of nitrate in most of the water, but can determine the concentration (mg / l N) of nitrite and then deduct the value of the content (in mg / l N) nitrate.

Please note that the refill kit's reference is 1MT333.

I-303 - Nitrates : 0.6 - 23.0 mg/L NO₃⁻-N
 METHOD COMPATIBLE WITH SEA WATER

Photopod
LS/SP

Reagent kit: 1MT184
 Preparation time: ~ 10 min

REAGENTS

Nitrates 1 Reagent	1RN015
Powder nitrate Reagent	1PN010
Demineralized water	1ED020

EQUIPMENT

Graduated Plastic Tube	14TP00 x 2
Crushing Rod	1AP018
Glass Tube	1CR099
Syringe 10ml	1SU013
Syringe 1 ml	1SU010
Filter holder	14PF09
Filter Paper	14PF05
Clamp	1PM010 (<i>consult us</i>)
Plastic spoon	1J000

TEST INSTRUCTIONS

Take the filter holder, unscrew it and introduce the filter paper, screw the bracket.
 In the first graduated plastic tube take, with the syringe introduce a 1 ml sample of water to analyse and complete with Demineralized water up to 10 ml
 Shake
 Add 1 spoonfull of Powder nitrate Reagent
 Close and shake 1 minute
 With the syringe 10ml take totally the sample (water+zinc)
 Fit the 10 ml syringe containing the sample on the filter holder and gently squeeze out a few drops of (=rinsing the filter).
 Filter the sample and introduce it in the second plastic tube up to 5 ml graduation.
 Add nitrate 1 Reagent up to 10 ml and shake.
 Wait 3 minutes.
 Proceed to the measurement
 Unscrew the filter holder, remove the soiled filter and clean all the water.

MEASUREMENT

Select the analysis **303 NO₃-N: 0.6 - 23.0 mg/L** (*Result in mg/L of N*)
 In the graduated plastic tube introduce 1 ml sample of water to analyze
 Fill the Tube to the 10 ml mark with Demineralized water.
 Fill a glass tube with this preparation using the plastic funnel then cover the tube.
 Put the black cover on top of the tube and press the key « zero »
 Remove the tube and put the sample tube to analyze
 Put the black cover on top of the tube and press the key « measure »

Nitrite correction

The method optionally reacts with the nitrite present in the sample. Nitrite content is low compared to that of nitrate in most of the water, but can determine the concentration (mg / l N) of nitrite and then deduct the value of the content (in mg / l N) nitrate.

Please note that the refill kit's reference is 1MT350.

I-305 - Nitrates : 0,06 - 1,80 mg/L NO₃⁻-N
 METHOD COMPATIBLE WITH SEA WATER



Reagent kit: 1MT184
 Preparation time: ~ 10 min

REAGENTS

Nitrates 1 Reagent	1RN015
Powder nitrate Reagent	1PN010

EQUIPMENT

Graduated Plastic Tube	14TP00 x 2
Crushing Rod	1AP018
Glass Tube	1CR099
Syringe 10ml	1SU013
Filter holder	14PF09
Filter paper	14PF05
Clamp	1PM010 (<i>consult us</i>)
Plastic spoon	1J0000

TEST INSTRUCTIONS

Take the filter holder, unscrew it and introduce it with the clamp, screw the bracket.
 In the first graduated plastic tube, take a 10 ml sample of water to analyze
 Add 1 spoonfull of Powder nitrate Reagent
 Close and shake 1 minute
 With the syringe 10ml take totaly sample (water+zinc)
 Fit the 10 ml syringe containing the sample on the filter holder and gently squeeze out a few drops of (=rinsing the filter).
 Filter the sample and introduce it in the second plastic tube graduated to 5 ml graduation.
 Add nitrate 1 Reagent to 10 ml and shake.
 Wait 3 minutes.
 Proceed to the measurement
 Unscrew the filter holder, remove the soiled filter and clean all the water.

MEASUREMENT

Select the analysis **305 NO₃-N: 0,06 - 1.80 mg/L** (*Result in mg/L of N*)
 Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
 Put the black cover on top of the tube and press the key « zero »
 Remove the tube and put the sample tube to analyze
 Put the black cover on top of the tube and press the key « measure »

Concentration as mg/L NO₃⁻ = result x 4,4

Nitrite correction

The method optionally reacts with the nitrite present in the sample. Nitrite content is low compared to that of nitrate in most of the water, but can determine the concentration (mg / l N) of nitrite and then deduct the value of the content (in mg / l N) nitrate.

Please note that the refill kit's reference is 1MT350.

I-306 - Nitrate : 0,1 - 20 mg/L NO₃⁻-N

Reagents kit reference: FTI2420702
Preparation time: ~ 5min

Photopod
LS/SP

REAGENTS

Reaction tube
Blank Tube (red label)
Nitrate-111

RECOMMENDED EQUIPMENT (consult us)

Automatic pipette 0,1 - 1 ml	1PA022
Pipette Tip 0,1 - 1 ml	1EU012
24 tubes stand Ø16	1PT013

TEST INSTRUCTIONS

With the pipette, put 0,5 ml of water to analyze in the reaction tube, close and gently shake.
Be careful, the reaction is exothermic and the tube becomes hot.
Add 0,2 ml of nitrate 111, close and return several times the tube.
Wait 5 minutes.

MEASUREMENT

Select the analysis **306 NO₃ -N : 0,02 - 20 mg/L**
Take the tube for the blank (tube with red label) and insert it in the photometer.
Put the black cover on top of the tube and press the key « zero ».
Remove the tube and put the sample tube to analyze.
Put the black cover on top of the tube and press the key « measure ».

Concentration as mg/L NO₃⁻ = result x 4,4

INTERFERENCES

Nitrite concentrations greater than 2 mg/L NO₂⁻ lead to higher test results. Great quantities of COD lead to higher test results.

I-320 - Nitrites: 0.01 - 0.60 mg/L NO₂-N
METHOD COMPATIBLE WITH SEA WATER

Reagent kit: 1MT027
Preparation time: ~ 6min



REAGENTS

Concentrated Ammoniac	1AC030
Z Indicator	1IZ000

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the glass tube
Add 7 drops of Z Indicator
Close the tube and shake.
Wait 5 minutes.
Add 7 drops of Concentrated Ammoniac
Close the tube and shake: a yellow color appears instantaneously
Proceed to the measurement

MEASUREMENT

Select the analysis **320 NO₂ -N: 0,01 - 0,60 mg/L** (Result in mg/L of N)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

$$\text{Concentration as mg/L NO}_2 = \text{result} \times 3,3$$

Please note that the refill kit's reference is 1MT334.

I-321- Nitrites: 0.01- 0.60 mg/L NO₂ -N

Reagent kit: 1MT165
Preparation time: ~ 11min

Photopod

SP

REAGENTS

Nitricol 1 Tablet 1AP109

EQUIPMENT

Graduated Plastic Tube 14TP00
Crushing Rod 1AP018
Glass Tube 1CR099
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 nitricol tablet crush it with the crushing rod and shake to dissolve ~ 1 min
Wait 10 min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **321 NO₂ -N: 0,01 - 0,60 mg/L** (*Result in mg/L of N*)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Concentration as mg/L NO₂ = result x 3,3

Please note that the refill kit's reference is 1MT335.

I-322- Nitrites: 0.4- 41 mg/L NO₂ -N

Reagent kit: 1MT166
Preparation time: ~ 3min

Photopod

SP

REAGENTS

Nitriphot 1 Tablet	1AP260
Nitriphot 2 Tablet	1AP260

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Nitriphot 1 tablet crush it with the crushing rod.
Close the tube and shake strongly to dissolve ~ 1-2 min
Add 1 Nitriphot 2 tablet crush it with the crushing rod and shake to dissolve ~ 15s
Wait 1 min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT


Select the analysis **322 NO₂ -N: 0.4 - 41 mg/L** (*Result in mg/L of N*)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Concentration as mg/L NO₂ = result x 3,3

Please note that the refill kit's reference is 1MT336.

I-323 - Nitrites: 4- 410 mg/L N

Reagent kit: 1MT166
Preparation time: ~ 3min

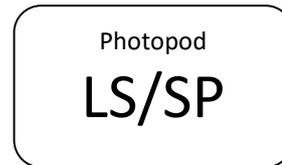

REAGENTS

Nitriphot 1 Tablet	1AP260
Nitriphot 2 Tablet	1AP260
Demineralized water	1ED010

EQUIPMENT

Graduated Plastic Tube	14TP00 x2
Crushing Rod	1AP018
Glass Tube	1CR099
Syringe 1 ml	1SU010
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

With the syringe, take a 1 ml sample of water to analyze, introduce it in the graduated tube then complete up to 10 ml with demineralized water.

Close and shake

Add 1 Nitriphot 1 tablet; crush it with the crushing rod.

Close the tube and shake strongly to dissolve ~ 1-2 min

Add 1 Nitriphot 2 tablet; crush it with the crushing rod and shake to dissolve ~ 15s

Wait 1 min

Fill a glass tube with this preparation using the plastic funnel then cover the tube.

Proceed to the measurement

MEASUREMENT

Select the analysis **323 NO₂ -N: 4 - 410 mg/L** (*Result in mg/L of N*)

In the graduated plastic tube introduce 1 ml sample of water to analyze

Fill the Tube to the 10 ml mark with Demineralized water.

Fill a glass tube with this preparation using the plastic funnel then cover the tube.

Put the black cover on top of the tube and press the key « zero »

Remove the tube and put the sample tube to analyze


Put the black cover on top of the tube and press the key « measure »

Concentration as mg/L NO₂ = result x 3,3

Please note that the refill kit's reference is 1MT336.

I-324 - Nitrite LR: 0,01- 1 mg/L NO₂⁻-N

Reagents kit reference: FTI2419018
Preparation time: ~ 12min

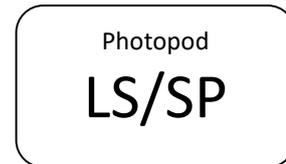

REAGENTS

Reaction tubes
Blank Tube (red label)
Nitrite-101

RECOMMENDED EQUIPMENT (consult us)

Automatic Pipette 1 - 5 ml	1PA023
Pipette Tip 1 - 5 ml	1EU013
24 tubes stand Ø16	1PT013

TEST INSTRUCTIONS

With the pipette, put 2 ml of water to analyze in the reaction tube, close and shake several times.

Add 1 level scoop of No. 8 (black) nitrite-101, close, and shake for 30 seconds.
Wait 5 minutes.

MEASUREMENT

Select the analysis **324 NO₂ -N: 0,01 - 1 mg/L** (Result in mg/L of N)
Take the tube for the blank (tube with red label) and insert it in the photometer.
Put the black cover on top of the tube and press the key « zero ».
Remove the tube and put the sample tube to analyze.
Put the black cover on top of the tube and press the key « measure ».

Concentration as mg/L NO₂ = result x 3,3

I-325 - Nitrite HR: 0,1- 5 mg/L NO₂⁻-N

Reagents kit reference: FTI2419018

Preparation time: ~ 12min

REAGENTS

Reaction tubes

Blank Tube (red label)

Nitrite-101

RECOMMENDED EQUIPMENT

Automatic pipette 0,1 - 1 ml	1PA022
Pipette Tip 0,1 - 1 ml	1EU012
24 tubes stand Ø16	1PT013

TEST INSTRUCTIONS

With the pipette, put 0,5 ml of water to analyze in the reaction tube, close and shake several times.

Add 1 level scoop of No. 8 (black) nitrite-101, close, and shake for 30 seconds.

Wait 5 minutes.

MEASUREMENT

Select the analysis **325 NO₂ -N: 0,1 - 5mg/L** (Result in mg/L of N)

Take the tube for the blank (tube with red label) and insert it in the photometer.

Put the black cover on top of the tube and press the key « zero ».

Remove the tube and put the sample tube to analyze.

Put the black cover on top of the tube and press the key « measure ».

Concentration as mg/L NO₂ = result x 3,3

I-331- Ozone: 0.30 - 4.00 mg/L O₃

Reagent kit: 1MT029
Preparation time: ~ 6min



REAGENTS

DPD 4 Tablet	1D4004P
Glycine Tablet	1NP000

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

1- TOTAL CHLORINE + OZONE

Take a 12,5 ml sample of water to analyze in the graduated plastic tube.
Add 1 DPD 4 tablet and shake to dissolve.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes.
Proceed to the measurement.
This gives the value 1: total chlorine + ozone in mg/l of O₃

2- TOTAL CHLORINE ONLY

Take a 12,5 ml sample of water to analyze in the graduated plastic tube.
Add 1 DPD 4 tablet and shake to dissolve.
Add 1 DPD Glycine tablet and shake to dissolve.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes.
Proceed to the measurement.
This gives the value 2: total chlorine in mg/l of O₃

3- OZONE

Concentration in mg/l of O₃ = value 1 - value 2

MEASUREMENT

Select the analysis **331 O₃:0.30 - 4.00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure ».

Please note that the refill kit's reference is 1MT337.

I-332 -Ozone: 0.03 - 0.65 mg/L O₃

Reagent kit: 1MT029
Preparation time: ~ 4min



REAGENTS

DPD 4 Tablet	1D4004P
Glycine Tablet	1NP000

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

1- TOTAL CHLORINE + OZONE

Take a 20 ml sample of water to analyze in the graduated plastic tube.
Add 1 DPD 4 tablet and shake to dissolve.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement
This gives the value 1: total chlorine + ozone in mg/l of O₃

2- TOTAL CHLORINE ONLY

Take a 20 ml sample of water to analyze in the graduated plastic tube.
Add 1 DPD 4 tablet and shake to dissolve.
Add 1 DPD Glycine tablet and shake to dissolve.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 2 minutes.
Proceed to the measurement
This gives the value 2: total chlorine in mg/l of O₃

3- OZONE

Concentration in mg/l of O₃ = value 1 - value 2


MEASUREMENT

Select the analysis **332 O₃:0.03 - 0.65mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT337.

I-340-Hydrogen Peroxide: 2 - 200 mg/l H₂O₂

Reagent kit: 1MT148
Preparation time: ~ 1.5min

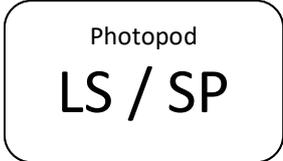

REAGENTS

Acidifying PT Tablet	1AP105
Hydrogen Peroxyde HR Tablet	1AP105

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Acidifying PT tablet; crush it with the crushing rod and shake to dissolve ~ 30s
Add 1 Hydrogen Peroxyde HR tablet; crush it with the crushing rod and shake to dissolve ~ 30s
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **340 H₂O₂ : 2 - 200mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT321.

I-341-Hydrogen Peroxyde: 0.05 - 2.00 mg/l H₂O₂

Reagent kit: 1MT149

Preparation time: ~ 2.5min

REAGENTS

Hydrogen Peroxyde LR Tablet 1AP104

EQUIPMENT

Graduated Plastic Tube 14TP00
Crushing Rod 1AP018
Glass Tube 1CR099
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Hydrogen Peroxyde LR tablet; crush it with the crushing rod and shake to dissolve ~ 30s
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 2 minutes
Proceed to the measurement

MEASUREMENT

Select the analysis **341 H₂O₂ : 0.05 - 2.00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT322.

I-350 - pH 6.8 - 8.6

Reagent kit: 1MT036
Preparation time: ~ 2min

Photopod

LS

REAGENTS

Phenol red 1PR008

EQUIPMENT

Glass Tube 1CR099
Syringe 10 ml 1SU013

TEST INSTRUCTIONS

With the syringe, take a 10 ml sample of water to analyze in the glass tube.
Add 16 drops of Phenol red.
Close and shake.
Proceed to the measurement

MEASUREMENT

Select the analysis **350 pH: 6.8 - 8.6mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT338.

I-392 - Phosphates: 0,50 - 13,0 mg/L PO₄³⁻P

I-380 - Phosphates: 1.00- 36.00 mg/L P₂O₅

Reagent kit: 1MT030

Preparation time: ~ 12min

Photopod

LS

REAGENTS

Phosphate 1 Reagent	1RP018
Phosphate 2 Reagent	1RP019

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

With the syringe, take a 10 ml of water to analyze in a glass tube

Add 8 drops of Phosphate 1 Reagent

Close and shake.

Add 8 drops of Phosphate 2 Reagent

Close and shake.

Wait 10 minutes.

Proceed to the measurement

MEASUREMENT

Select the analysis **392 PO₄ -P: 0.50 - 13.0 mg/L** (Result in mg/L of P)

or **380 P₂O₅ : 1.00- 36.0 mg/L** (Result in mg/L of P₂O₅)

Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer

Put the black cover on top of the tube and press the key « zero »

Remove the tube and put the sample tube to analyze

Put the black cover on top of the tube and press the key « measure »

Concentration as mg/L of PO₄ = result as mg/L of P x 3,1

Please note that the refill kit's reference is 1MT352.

I-390 - Phosphates: 0.06 - 1.30 mg/L PO₄³⁻-P

Reagent kit: 1MT186
Preparation time: ~ 5 min

Photopod

SP

REAGENTS

Phosphate 1 Tablet	1AP177
Phosphate 2 Tablet	1AP177

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyse in the graduated plastic tube.
Add 1 Phosphate 1 Tablet, crush it with the crushing rod.
Close the tube and shake to dissolve ~ 3 min
Add 1 Phosphate 2 Tablet, crush it with the crushing rod and shake to dissolve ~ 1 min
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 1 minute.
Proceed to the measurement

MEASUREMENT

Select the analysis **390 PO₄-P 0.65 - 1.30 mg/L** (*Result in mg/L of P*)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Concentration as mg/L of PO₄ = result as mg/L of P x 3,1

Please note that the refill kit's reference is 1MT354.

I-391 - Phosphates: 0.06 - 1.60 mg/L PO₄³⁻-P

Reagent kit: 1MT030
Preparation time: ~ 12min

Photopod

LS

REAGENTS

Phosphate 1 Reagent	1RP018
Phosphate 2 Reagent	1RP019

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

With the syringe, take a 10 ml of water to analyze in a glass tube
Add 2 drops of Phosphate 1 Reagent
Close and shake.
Add 2 drops of Phosphate 2 Reagent
Close and shake.
Wait 10 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **391 PO₄-P: 0.06- 1.60 mg/L** (Result in mg/L of P)
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Concentration as mg/L of PO₄ = result as mg/L of P x 3,1

Please note that the refill kit's reference is 1MT352.

I-402 - Total Phosphates: 0,05 - 3 mg/L PO₄-P

Reagents kit reference: FTI2419019
Preparation time: ~ 40 min

Photopod
LS / SP

REAGENTS

Reaction Tubes
Blank Tube (red label)
Phosphate-101
Phosphate-102
Phosphate-103

RECOMMENDED EQUIPMENT (consult us)

Automatic Pipette 1 - 5 ml	1PA023
Pipette Tip 1 - 5 ml	1EU013
24 tubes stand Ø16	1PT013
Wooden clamp	1PT007
Heating reactor	1RD010

TEST INSTRUCTIONS

Turn on the heating reactor. Preheat at 100 °C.
Take one reaction tube and with the pipette, put 5 ml of water to analyse.
Add 1 level scoop of No. 4 (white) phosphate-103, close immediately, and shake for 30 seconds.
Put the tubes in the heating reactor at 100°C during 30 minutes.
After 30 minutes, take the tubes with the wooden clamp (be careful, they are very hot) and shake gently. Put the tube in the tube stand et let cool to room temperature (>20 minutes).
Add 2 drops (0,1 ml) phosphate-101, close and shake several times.
Add 1 level scoop of No. 4 (white) phosphate-102, close, and shake for 30 seconds.
Wait 5 minutes.

MEASUREMENT


Select the analysis **402 PO₄-P: 0,05 - 5mg/L** (Result in mg/L of P)
Take the tube for the blank (tube with red label) and insert it in the photometer.
Put the black cover on top of the tube and press the key « zero ».
Remove the tube and put the sample tube to analyze.
Put the black cover on top of the tube and press the key « measure ».

Concentration as mg/L of PO₄ = result as mg/L of P x 3,1

NOTA - If the analysis is performed without digestion only PO₄ ions are determined.

I-410 - Potassium : 2.00 - 15.0 mg/L K

Reagent kit: 1MT168
Preparation time: ~ 4min

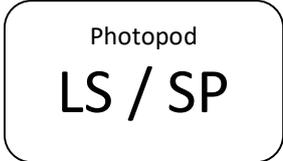

REAGENTS

Potassium Tablet 1AP189

EQUIPMENT

Graduated Plastic Tube 14TP00
Crushing Rod 1AP018
Glass Tube 1CR099
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Potassium tablet, crush it with the crushing rod and shake to dissolve ~ 45s
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 3 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **410 K: 2.00- 15.0mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT340.

I-420 - Silica: 10 - 300 mg/L SiO₂
I-421 - Silica: 0,20 - 10 mg/L SiO₂

Photopod

LS

Reagent kit: 1MT040
Preparation time: ~8min

REAGENTS

Ammonium Molybdate	1AM010
Sulfuric Acid ¼	1AS013
Oxalic Acid 10%	1AO000

EQUIPMENT

Graduated Plastic Tube	14TP00
Plastic spoon	1J0000
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 20 ml sample of water to analyze in the graduated plastic tube.
Add 1 spoonful to the brim of Ammonium Molybdate and shake to dissolve 30 s
Add 7 drops of Sulphuric acid ¼ and shake.
Wait 5 minutes.
Add 15 drops of Oxalic acid 10% and shake.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **420 SiO₂: 10 - 300mg/L**
or **421 SiO₂: 0.20 - 10 mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT341.

I-422 - Silica: 5 - 150 mg/L SiO₂

Reagent kit: 1MT173
Preparation time: ~ 12 min

Photopod

SP

REAGENTS

Silica HR 1 Tablet	1AP290
Silica PR Tablet	1AP290
Silica HR 2 Tablet	1AP290

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Silica HR 1 Tablet, crush it with the crushing rod. Close the tube and shake strongly to dissolve. ~ 2 min
Add 1 Silica HR 2 Tablet, crush it with the crushing rod. Close the tube and shake strongly to dissolve. ~ 2 min
Add 1 Silica PR Tablet; crush it with the crushing rod. Close the tube and shake strongly to dissolve ~ 6 min
Wait 2 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **422 SiO₂ : 5 - 150mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT342.

I-423 - Silica: 0.05 - 10 mg/L SiO₂

Reagent kit: 1MT170
Preparation time: ~ 12 min



REAGENTS

Silica 1 Tablet	1AP181
Silica PR Tablet	1AP181
Silica 2 Tablet	1AP181

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Silica 1 Tablet; crush it with the crushing rod. Close the tube and shake strongly to dissolve. ~ 2 min
Add 1 Silica PR Tablet; crush it with the crushing rod. Close the tube and shake strongly to dissolve. ~ 6 min
Add 1 Silica 2 Tablet; crush it with the crushing rod. Close the tube and shake strongly to dissolve. ~ 2 min
Wait 2 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement

MEASUREMENT

Select the analysis **423 SiO₂ : 0.05 -10mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT343.

I-430 - Sulfates: 10 - 400 mg/L SO₄²⁻

Reagent kit: 1MT080
Preparation time: ~ 11min

Photopod

LS

REAGENTS

Sulfates n°1 Reagent	1RS015
Sulfates n°2 Reagent	1RS016

EQUIPMENT

Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

With the syringe, take a 10 ml sample of water to analyze in the glass tube
Add 5 drops of Sulfates n°1 Reagent, close the tube and shake strongly to 15s
Add 10 drops of Sulfates n°2 Reagent, close the tube and shake strongly to 15s.
Wait 10 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **430 SO₄: 10 - 400mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT344.

I-431-Sulfates: 10-200 mg/L SO₄²⁻

Reagent kit: 1MT171
Preparation time: ~ 6min

Photopod

SP

REAGENTS

Sulfate Tablet 1AP154

EQUIPMENT

Graduated Plastic Tube 14TP00
Crushing Rod 1AP018
Glass Tube 1CR099
Plastic funnel 1EP021 (*consult us*)

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Sulfate Tablet; crush it with the crushing rod and shake to dissolve ~ 45s
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Wait 5 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **431 SO₄: 10 - 200mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT041.

I-440 - Sulfide: 0.05 - 0.60 mg/L S

Reagent kit: 1MT172
Preparation time: ~ 6 min

Photopod

LS

REAGENTS

Sulfide n°1 Tablet	1AP168
Sulfide n°2 Tablet	1AP168

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099

TEST INSTRUCTIONS

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Sulfide n°1 Tablet and 1 Sulfide n°2 Tablet; Crush it with the crushing rod and shake to dissolve. ~ 1 min
Wait 5 minutes.
Proceed to the measurement

MEASUREMENT

Select the analysis **440 S: 0.05- 0.60mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT345.

I-460 - Zinc: 0.05- 4.00 mg/l Zn

Reagent kit: 1MT190
Preparation time: ~ 2min

Photopod

LS

REAGENTS

Zinc 1 Reagent	1RZ011
Zinc 2 Reagent	1RZ012

EQUIPMENT

Crushing Rod	1AP018
Glass Tube	1CR099
Syringe 10 ml	1SU013

TEST INSTRUCTIONS

With the syringe, take a 10 ml sample of water to analyze and introduce in the glass tube
Add 5 drops of Zinc 1Reagent
Shake
Wait 1 minute
Add 10 drops of Zinc 2 Reagent
Shake
Proceed to the measurement

MEASUREMENT

Select the analysis **460 Zn: 0.05- 4.00mg/L**
Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer
Put the black cover on top of the tube and press the key « zero »
Remove the tube and put the sample tube to analyze
Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT356.

I-461 - Zinc: 0.10- 4.00 mg/IZn

Reagent kit: 1MT043
Preparation time: ~ 6 min



REAGENTS

Pack Zinc	1PZ001
Zinc Tablets	1PZ001
Dechlor Tablets	1PZ001
EDTA Tablets	1PZ001

EQUIPMENT

Graduated Plastic Tube	14TP00
Crushing Rod	1AP018
Glass Tube	1CR099
Plastic funnel	1EP021 (<i>consult us</i>)

TEST INSTRUCTIONS

Water not containing copper or chlorine

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Zinc tablet, crush it with the crushing rod and shake to dissolve.
Wait 5 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement.

Water containing copper

Follow the instructions for « Water not containing copper or chlorine» above and proceed to the measurement.
The result is the concentration in zinc and copper $Conc(Zn + Cu)$
Transfer the content of the glass tube in the graduated plastic tube.
Add 1 EDTA tablet, Add 1 Zinc tablet, crush it with the crushing rod and shake to dissolve. (the color due to zinc disappear, the color due to copper remains).
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement
The result is the concentration in copper $Conc(Cu)$

Concentration zinc is :

$$Conc(Zn) = Conc(Zn + Cu) - Conc(Cu)$$

Water containing chlorine

Take a 10 ml sample of water to analyze in the graduated plastic tube.
Add 1 Dechlor tablet, crush it with the crushing rod and shake to dissolve.
Add 1 Zinc tablet, crush it with the crushing rod and shake to dissolve.
Wait 5 minutes.
Fill a glass tube with this preparation using the plastic funnel then cover the tube.
Proceed to the measurement.

MEASUREMENT

Select the analysis **461 Zn: 0,10 - 4,00 mg/L**

Fill a glass tube with water to analyze without reagent, put the cap and insert it in the photometer.

Put the black cover on top of the tube and press the key « zero ».

Remove the tube and put the sample tube to analyze.

Put the black cover on top of the tube and press the key « measure »

Please note that the refill kit's reference is 1MT346.