

Aquinity² P10

Produces Ultra Pure Water




WATER
PURIFIER






WATER
PURIFIER

The Aquinity² P10 system is designed for the production of ultra pure water (0.055 μ S/cm) out of tap water.

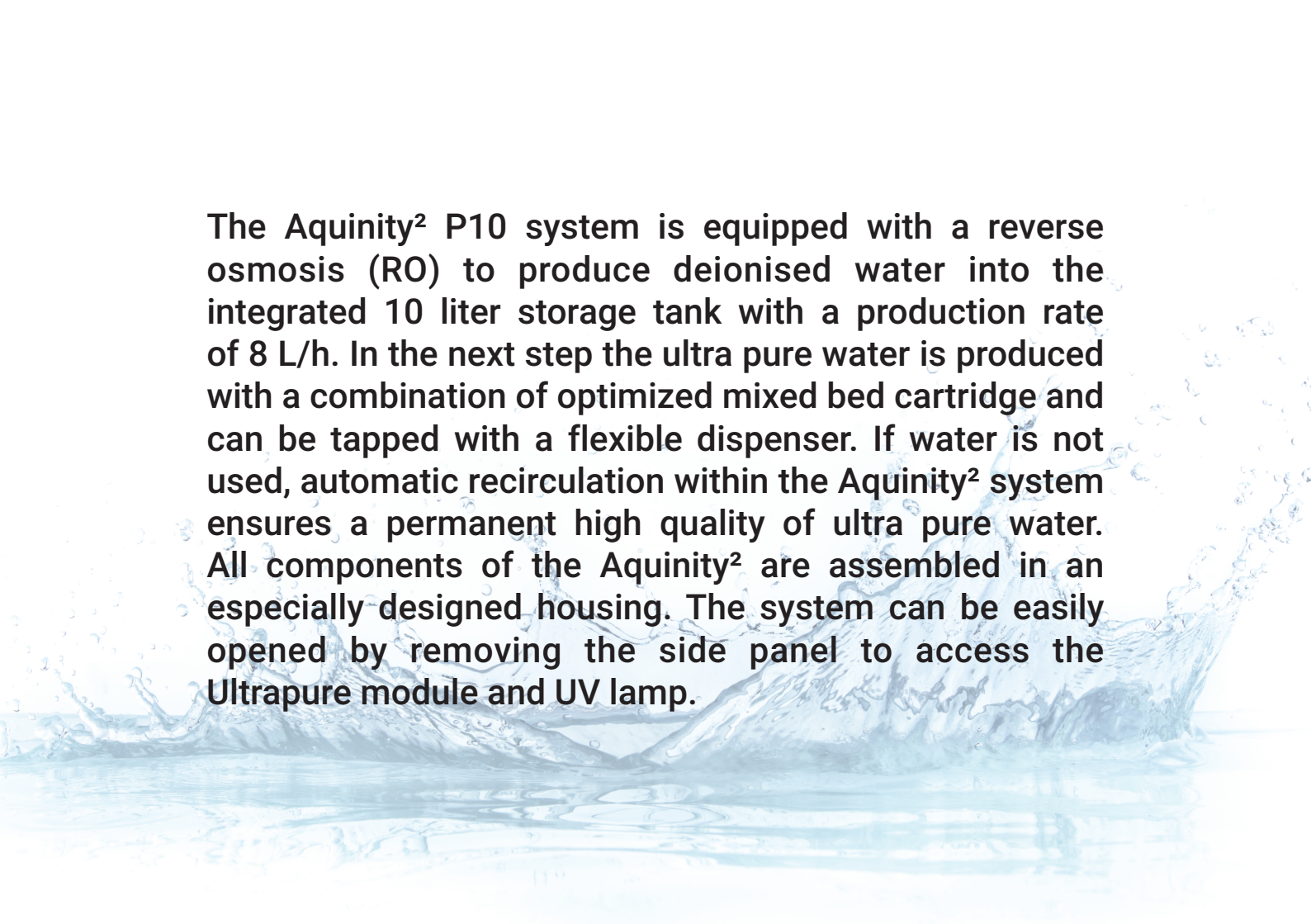A dynamic splash of water with many droplets and ripples, creating a sense of movement and freshness. The water is light blue and white, set against a plain white background.

The Aquinity² P10 system is equipped with a reverse osmosis (RO) to produce deionised water into the integrated 10 liter storage tank with a production rate of 8 L/h. In the next step the ultra pure water is produced with a combination of optimized mixed bed cartridge and can be tapped with a flexible dispenser. If water is not used, automatic recirculation within the Aquinity² system ensures a permanent high quality of ultra pure water. All components of the Aquinity² are assembled in an especially designed housing. The system can be easily opened by removing the side panel to access the Ultrapure module and UV lamp.

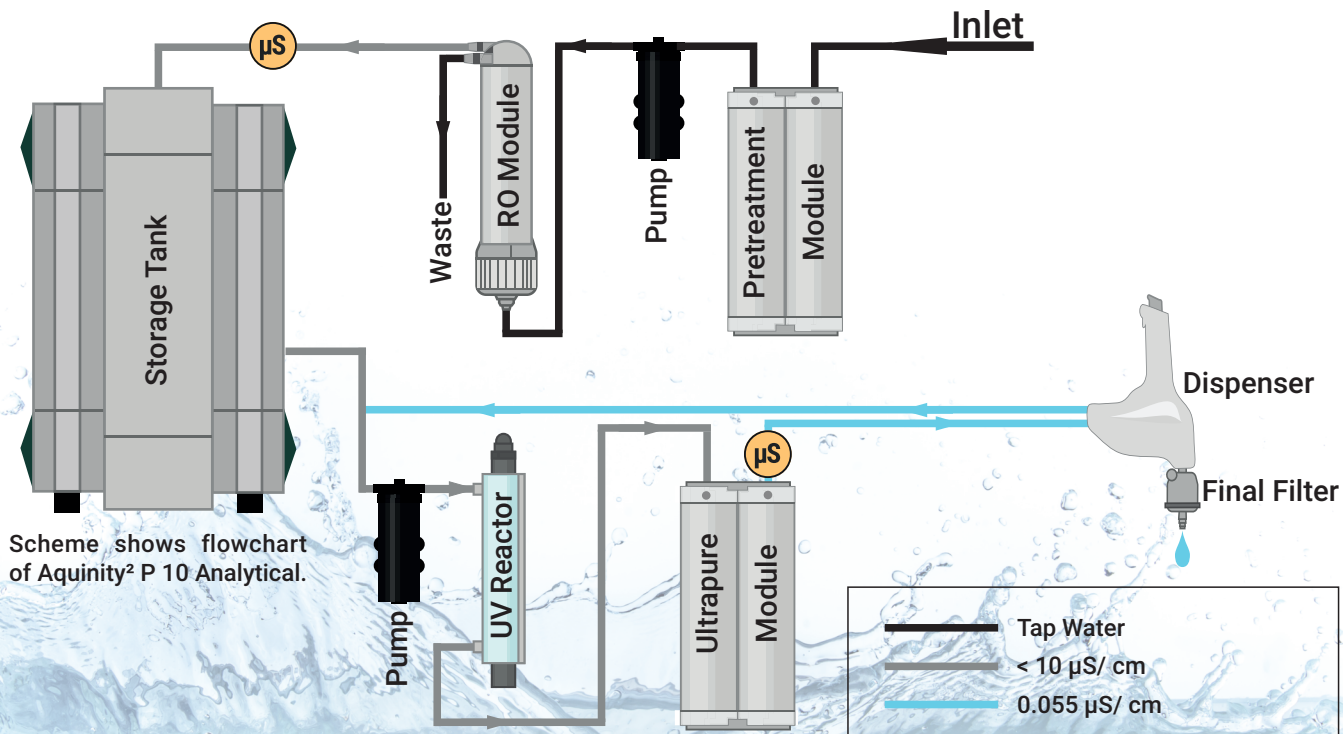

Technical Specifications



ultra pure water quality	0.055 $\mu\text{S}/\text{cm}$; Type I
resistivity	18.2 $\text{M}\Omega\cdot\text{cm}$
total organic carbon (TOC)	< 10 ppb Reagent < 3 ppb Analytical
dispensing flow rate	1.5 L/min
productivity rate	8 L/h*
bacteria	< 1 cfu/mL**
particulate	> 0.2 μm less than 1 particulate/mL
RNAse	< 1 pg/mL**
DNAse	< 5 pg/mL**
dimensions, weight, power	504 x 340 x 535 mm, 16-20 kg, 110-230 V

* not possible to tap pure water ** with microbiological final filter

Feed Water Requirements

type of feed water	tap water
feed water conductivity	< 1400 $\mu\text{S}/\text{cm}$
inlet pressure	1.5 to 6 bar
free chlorine	< 0.1 mg/L
Silt Density Index (SDI)	< 3
pH	3 to 9
temperature	5 to 25 °C

Configurations



Aquinity²P10 is available in two configurations (Reagent & Analytical) to fit the specific requirements of ultra pure water quality for different applications.

model	UV reactor	Cat. No.
Reagent	-	114-0070
Analytical	+	114-0071

Options

µS-control

The µS-control checks the conductivity of feed water to protect the cartridges against inappropriate feed water. If conductivity is too high the water flow will be rejected.

TOC monitoring (TI version)

The TOC monitoring during production and intermittent measurements during non-use periods allows to check the organic content in water continuously between 1 and 999 ppb.

volumetric dispensing

Our dispenser allows the volumetric controlled dispensing of water with an increment of 0.1 L and a tap volume from 0.1 to 99 L. The system prevents overflow of containers and allows to dispense water without supervising.




bench integrated (BI version)

To save valuable benchtop space the Aquinity² can be ordered as BI configuration. The system will be installed underneath the work bench and only the display and dispenser are mounted on the wall.




User Interface & Software


Currently recorded data and warning messages will be displayed on the touch screen monitor. The software furthermore supports the user with a self-diagnostic module which reduced service time and costs.

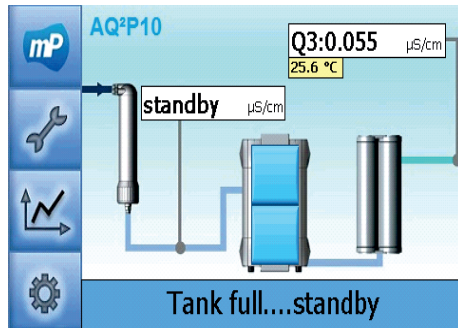


Main Screen: The Software allows the user to see all information to use maintenance tool and track back historical values up to 1 year.

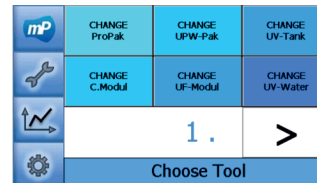

Helping tools to maintain system

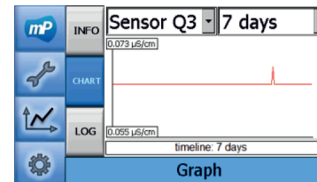

Trackable values up to 1 year

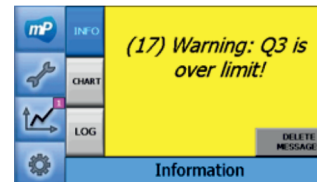

Warning of customer in case of over limit and service case



Consumables

Description	Cat. No.
Ultrapure Module MemPak LS (Life Science)	190-0087
Ultrapure Module MemPak AL (Analytical)	190-0088
UV lamp	921-0508
Pretreatment Module ProPak R10	290-0065
Final Filter, capsule 0.2 µm	190-0013*
Disinfection Tablets	290-0227

* 190-0082 Final Filter to reduce Edotoxin, DNA + RNase



WATER
PURIFIER

Made in Germany, more than 25 years...



membraPure GmbH
Wolfgang-Küntscher-Str. 14,
16761 Hennigsdorf/Berlin, Germany
Tel.: +49 3302 - 201 20 0



info@membrapure.de
www.membrapure.de